



Experience visible lift,
natural rejuvenation, and
timeless confidence.

www.aestheon.com/ullift

Ullift
Redefining the Art of Lifting

Why Ullifting?

Ullift® introduces an innovative dual-mode system that combines Linear and Dot Modes in a single session—targeting both deep and superficial dermal layers simultaneously. By stimulating multilayer collagen regeneration, Ullift® delivers precise energy to the skin for natural lifting and firming effects. The result is a faster, safer, and more comfortable treatment with visible results and minimal downtime.



1.5 mm



3.0 mm



4.5 mm



6.0 mm



9.0 mm



• Forehead / Brow

• Periorbital / Eye


• Chin / Jaw

• Neck

• Décolleté

The Precise Dual HIFU Modes

Ulift® is a next-generation HIFU lifting device designed with two distinct treatment modes – Linear Mode and Dot Mode – to deliver precise energy to multiple dermal layers for optimal lifting and rejuvenation outcomes. Linear Mode employs a 1.5 mm cartridge to deliver continuous energy lines across the superficial dermal layer in a safe and controlled manner.

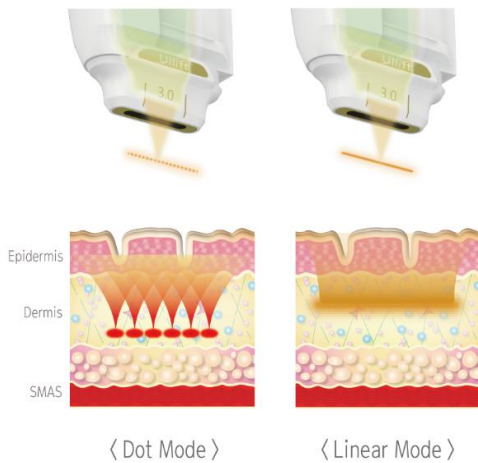

This technique enhances collagen neogenesis and fibrosis modulation, achieving greater efficacy compared to conventional HIFU systems while maintaining excellent patient comfort and minimal pain response. Dot Mode utilizes 3.0 mm, 4.5 mm, and 6.0 mm cartridges to deliver slightly higher-intensity energy to the deep dermal and SMAS layers. By stimulating collagen remodeling and fibrosis improvement within the deeper tissue structures, Dot Mode provides long-lasting lifting effects with mild discomfort and a high safety profile.

The Easiest GUI (Graphical User Interface)

Ulift® is equipped with a highly intuitive and simplified GUI, thoughtfully designed to ensure effortless operation for every practitioner. Clear icons, logical navigation, and a clean visual layout allow users to quickly understand the system and begin treatments with confidence. Even first-time operators can easily adjust parameters, select cartridges, and execute procedures without additional complexity or long learning curves. With Ulift®, usability becomes seamless—enabling more focus on patient care and treatment results rather than system control.



Premium & Emotional

With its expansive 12.1-inch screen, Ullift® enhances every detail for a clearer, more confident treatment experience.

Ullift® is designed to deliver both precision and comfort through its advanced user interface. Unlike conventional HIFU devices with smaller displays, Ullift® features a large, high-resolution LCD screen that enhances visibility and enables practitioners to monitor every detail with confidence. Its clean, intuitive graphical user interface (GUI) has been meticulously designed to support accurate operation and reduce the possibility of user error, even during fast-paced clinical sessions. The streamlined single-page layout presents all key parameters in one view, allowing practitioners to make real-time adjustments without unnecessary navigation. This efficient interface design not only ensures faster workflow and improved safety but also provides an intuitive experience where even first-time users can operate the system with ease — without the need for additional training.

Less Pain, Greater Comfort

Ullift® was developed with a strong focus on minimizing pain and maximizing patient comfort.

Traditional HIFU procedures often caused significant discomfort for both practitioners and patients. To address this, the Ullift® engineering team dedicated extensive research to reducing pain while enhancing treatment efficacy. Through its ultra-fast scanning speed of up to 0.5 second per shot and smooth, precisely controlled HIFU targeting, Ullift® successfully alleviates the severe pain often associated with conventional HIFU treatments. This breakthrough technology allows for faster procedures, gentler energy delivery, and greater satisfaction for both patients and practitioners alike.

Treatment Indications

- Lifting / tightening
- Face /jawline contouring
- Deep wrinkles and fine wrinkles
- Lifting around the eyes / eyelids
- Body contouring



Visible Clinical Improvements after Lifting, Tightening & Contouring

PRACTITIONER BENEFITS

- **Exceptional ROI** with proven reliability and long-term value
- **Bridges the gap** between traditional HIFU and next-generation lifting solutions
- **12.1-inch-wide display** and intuitive GUI for faster, safer operation
- **Efficient single-page** interface minimizing user error and improving workflow
- Designed for **consistent results and customer satisfaction**


PATIENT BENEFITS

- **Comfortable treatments** with reduced pain through smooth HIFU targeting
- **Visible lifting results** after a single session with no downtime
- **Multidepth energy delivery** stimulating collagen at both superficial and deep layers
- **Natural and long-lasting outcomes** with improved skin firmness and elasticity
- A **gentler, smarter HIFU experience** for every skin type

Jawline contouring



Eye wrinkles

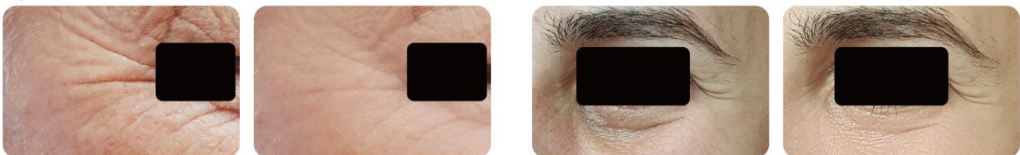

Skin tightening



Face lifting





Optimized
to enhance treatment efficacy
and patient comfort.



System Specifications

Type	HIFU(High Intensity Focused Ultrasound)
Handpiece	Line-type dual handpieces
Frequency	4 Mhz, 7 Mhz
Cartridge Depth	Face : 1.5, 3.0, 4.5 mm Body : 6.0, 9.0 mm
Treatment Length	25mm Max
Energy	Face 1.5J Max (0.1J per a step) Body 2.5J Max (0.1J per a step)
Pitch	1.0~2.0 mm (0.1 mm)
Interface	12.1 inch 1280*800 TFT LCD
Dimensions (WxDxH)	374 x 283 x 300(mm)
Input Power	AC 100 ~ 240V, 50/60Hz



Aestheon Co., Ltd.

1F, Centumdong-ro 45, Haeundae-gu, Busan 48059, Republic of Korea.

Tel: +82-70-4090-2027 / Fax: +82-51-742-7990

Website: www.aestheon.com